



MPR Pharmaceutical Refrigerators

+2°C to +23°C

684 L / 671 L (R type)

The MPR-722 (shelves) and MPR-722R (drawers) pharmaceutical refrigerators deliver stable, reliable temperature control for clinical, pharmaceutical, biomedical storage and processing. Fast recovery following door openings assures safe storage even with frequent access and high ambient conditions and to remove heat loads generated by powered instruments.

Airflow Assures Quick Recovery, Interior Uniformity

Vertical, forced air circulation with blower creates a uniform temperature at all shelf or drawer levels, top to bottom and front to back. Internal shelf or drawer clearances permit positive airflow to help restore temperature to setpoint following door openings, and safely accommodates busy, high traffic demand for access to stored contents in pharmacies, clinics and hospitals.

Precision Temperature Control Protects Inventory

Because many pharmaceuticals are degraded or destroyed if accidentally frozen, it is important that temperature be adjusted to desired setpoint in accordance with drug manufacturers' insert specifications. MPR Series pharmaceutical refrigerators operate over a temperature range of 2°C to 23°C. Actual interior air temperature is displayable in 0.1°C increments.

Choose Your Inventory Management Preference

The MPR Series pharmaceutical refrigerator offers two choices in inventory management. The MPR-722 includes 4 open-wire shelves strong enough to hold bulk loads. The MPR-722R is equipped with heavy-duty, fully-extendable painted steel pull-out drawers for more convenient access to the stored product such as large bottles or reagent kits.




Airflow Prevents Temperature Stratification

The top-mounted air circulation blower directs air downward so it goes along the back, side and front of the load-in drawers or shelves. Regardless of location, all products are subjected to the same safe temperature to assure repeatability of stored product viability.



Microprocessor Controls Simplify Operation

A microprocessor controller monitors all functions. Setpoint is factory set at 5°C. A soft key interface allows setpoint adjustment throughout the operating range. Audible and visual warnings indicate door ajar or temperature deviation from setpoint. The large easy-to-read LED display confirms actual temperature.

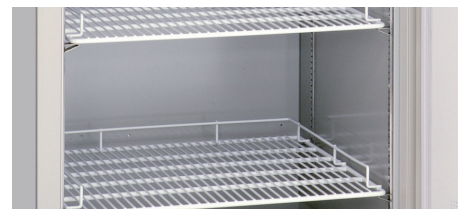

Adjustable Shelves Offer Flexibility in Storage Options

Open wire shelves (shown) are adjustable to allow for loads of varying heights such as bulk storage or process equipment such as fraction collectors. The Model MPR-722R with pull-out drawers (not shown) offers a practical solution with additional ergonomic advantages.

