

# Steam Sterilizer (Laboratory Use, Large Capacity)

## SQ500/510

Operating temp. range 45°C to 135°C

Internal capacity 47.5L

**High-performance, easy operation and ergonomically designed sterilizer featuring easy top loading and handling of samples**

- Maximum sterilizing temperature goes up to 135°C suitable for protein modification
- By simple setting and operation, these sterilizers manage ordinary sterilization, sterilization of culture mediums and liquids, and melting of culture mediums
- Custom programs can be saved for repeated use
- Standard equipped with slow release valve and cooling fan to shorten cool down time
- Mobile on wheels



47.5L  
SQ500

### Specifications

Product code	241008	241057
<b>Model</b>	<b>SQ500</b>	<b>SQ510</b>
System	Automatic high-pressure steam sterilization*1	
Operating temperature range	45~135°C (45~80°C (Preheating temperature)/45~60°C (Heat-retention process)/65~100°C (Solution process)/105~135°C (Sterilization) )	
Max. working pressure	0.255MPa	
Ambient temperature	5~35°C	
Lid	Manual upward opening with an interlock for safety	
Heater	1000Wx2 units	
Exhaust valve	One exhaust valve and one slow release valve	
Connection ports for optional accessories	Total 3 ports. Female thread for sample temp. sensor (1/4" ), Female thread for chamber temp. sensor (1/4" ), Female thread pressure sensor (Branching from the solenoid valve tubing)	
Cooling fan	Axial fan motor	
Temp. controller	PID control by microprocessor	
Temp. display / setting	Digital display / digital setting by ▲/▼ keys	
Timer / Timer resolution	0 or 1 min. to 99 hrs 59 min. / 1 minute	
Operation mode	Instrument sterilization course, Liquid sterilization course, Sterilization and keep warm course, Melting and keep warm course, Customer-programmed course	
Other functions	Key locking, Presetting, Saving, preheating, Forced cooling, Sample temperature sensor (option), Pattern locking, Up to 20 error log saving, Display of accumulated working time / Present time, ON-OFF beeping setting	
Safety device	Sensor failure detection, SSR short-circuit, Broken heater wire, Prevention of idle heating (Liquid expansion type), Alarm against the absence of a drain bottle, Failure in locking the lid, Memory error detection, Pressure relief valve (0.255 MPa)	
External dimensions*2	W520xD660xH846mm	
Internal dimensions of chamber	I.D.370xD442 mm	
Internal capacity	47.5L	
Weight	Approx. 95kg	
Power source (50/60Hz single phase)	AC100~120V (20.5~24.5A)	AC200~240V (10.5~12.5A)
Accessories	2 pcs. stainless steel mesh basket (Dia.332xD195.5mm) OSQ-90 Vapor cupx1, Drain bottlex1, Drain boardx1, Filterx1	


The length of the power cord is about 3m outside the unit.

\*1 This is not a medical device. \*2 Do not include protrusions


## Low Height Sterilizers



## Control panel

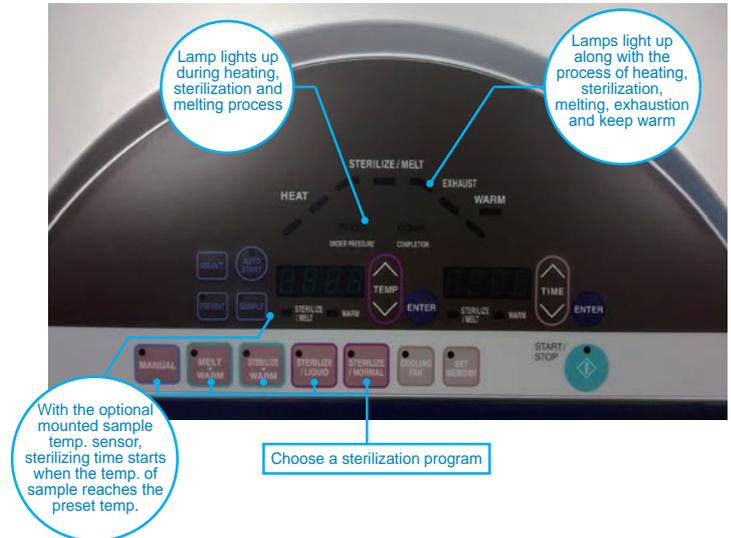

## Comparison of the cooling time of model SQ and conventional autoclave

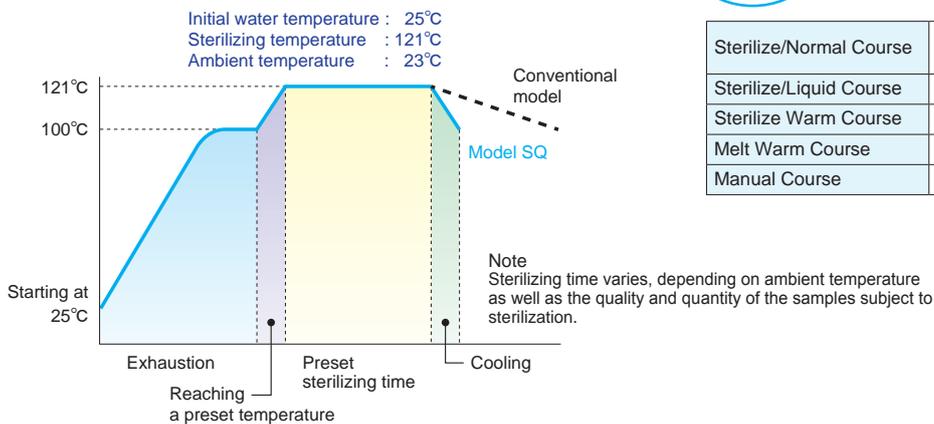

Sterilize/Normal Course	Sterilization of equipment such as flask, beaker, test tube, scissors
Sterilize/Liquid Course	Sterilization of culture and reagents and keep warm
Sterilize Warm Course	Dissolve and keep warm of the agar medium
Melt Warm Course	Sterilization of liquid, purified water and dilution water
Manual Course	Customized temperature and time settings

## Optional items

No.	Description	Model	Product code
①	Mesh basket	OSQ-90	241099
②	Mesh basket with 2 stacking fittings	OSR-40	241090
③	Mesh basket with 1 adjustable stainless steel perforated plate	OSR-50	241097
④	Stainless bucket	OSN14	241152
⑤	*Chamber temperature sensor, 3pcs/set	OSR-70	241153
⑥	*Sample temperature sensor, 1pc	OSN16	241154
⑦	*Pressure gauge	OSN18	241155
⑧	*External alarm output terminal	OSN20	241157
⑨	*Temperature output terminal (4~20mA)	OSN22	241158
⑩	*Time-up output terminal (Reray, contact output)	OSN24	241159

\* Please specify when ordering main unit.

## Stainless baskets



## Stainless buckets



# Steam Sterilizer (Laboratory Use, Large Capacity)

SQ500C/510C/810C



Operating temp. range 45°C to 135°C

Max. operational Pressure 0.255MPa

Internal capacity 50L SQ500C/510C 80L SQ810C

High-performance, easy operation and ergonomically designed sterilizer featuring easy top loading and handling of samples




- Maximum sterilizing temperature goes up to 135°C suitable for protein modification
- By simple setting and operation, these sterilizers manage ordinary sterilization, sterilization of culture mediums and liquids, and melting of culture mediums
- Custom programs can be saved for repeated use
- Standard equipped with timer start and pre-heating features for efficient use of time and cooling fan to shorten cool down time
- Mobile on wheels
- **Enhanced safety device**
  - Lid interlock mechanism
  - Drain bottle unset alarm
  - Memory malfunction
  - Automatic shutdown when malfunction
  - Self-diagnostic functions

## Specifications

Model	SQ500C	SQ510C	SQ810C
System	Automatic high-pressure steam sterilization		
Operating temperature range	45~135°C		
Max. working pressure	0.255MPa		
Ambient temperature	5~35°C		
Lid	Manual upward opening with an interlock for safety		
Heater	120V, 1440Wx2 units		115V, 1600Wx2 units
Exhaust valve	One exhaust valve and one slow release valve		
Connection ports for optional accessories	Total 3 ports. Female thread for sample temp. sensor (1/4" ), Female thread for chamber temp. sensor (1/4" ), Female thread pressure sensor (Branching from the solenoid valve tubing)		
Cooling fan	Axial fan motor		
Temp. controller	PID control by microprocessor		
Temp. display / setting	Digital display / digital setting by ▲/▼ keys		
Timer / Timer resolution	0 or 1 min. to 99 hrs 59 min. / 1 minute		
Operation mode	Instrument sterilization course, Liquid sterilization course, Sterilization and keep warm course, Melting and keep warm course, Customer-programmed course		
Other functions	Key locking, Presetting, Saving, preheating, Forced cooling, Sample temperature sensor (option), Pattern locking, Up to 20 error log saving, Display of accumulated working time / Present time, ON-OFF beeping setting		
Safety device	Sensor failure detection, SSR short-circuit, Broken heater wire, Prevention of idle heating (Liquid expansion type), Alarm against the absence of a drain bottle, Failure in locking the lid, Memory error detection, Pressure relief valve (0.255 MPa)		
External dimensions (WxDxHmm)	520x660x846		520x660x1161
Internal dimensions of chamber	I.D.370xD470 mm		I.D.370xD750 mm
Internal capacity	50L		80L
Weight	Approx. 95kg		Approx. 105kg
Power source	AC100~120V (20.5~24.5A)	AC200~240V (10.5~12.5A)	AC200~240V (10.5~12.5A)
Accessories	2 pcs. stainless steel mesh basket (Dia.332xD195.5mm), OSQ-90		3 pcs. stainless steel mesh basket (Dia.332xD195.5mm), OSQ-90
	Vapor cupx1, Drain bottlex1, Drain boardx1, Chemical indicator 1 set (30 pieces), Filterx1		

## Low Height Sterilizers

SQ500C(low height type)

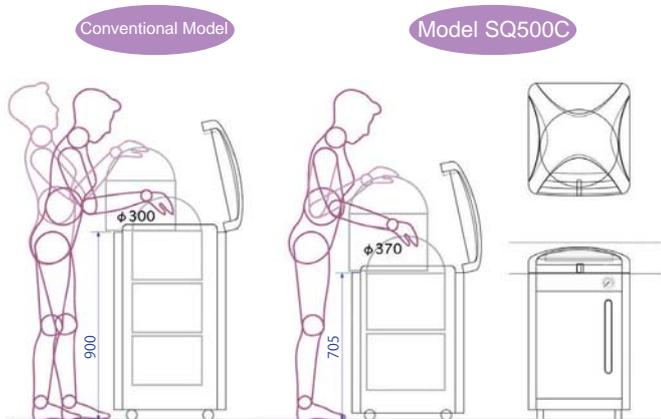

## Control panel



Sterilize/Normal Course	Sterilization of equipment such as flask, beaker, test tube, scissors
Sterilize/Liquid Course	Sterilization of culture and reagents and keep warm
Sterilize Warm Course	Dissolve and keep warm of the agar medium
Melt Warm Course	Sterilization of liquid, purified water and dilution water
Manual Course	Customized temperature and time settings

## Optional items

Baskets with adjustable stainless steel perforated plate

Mesh baskets

Stainless solid baskets

Baskets with stacking fittings



### Stainless baskets

### Stainless buckets

Product code	Model	Description
241099	OSQ-90	Mesh basket
241090	OSR-40	Mesh basket with 2 stacking fittings
241097	OSR-50	Mesh basket with 1 adjustable stainless steel perforated plate
241098	OSR-60 *Discontinued	Stainless solid basket
241152	OSN14	Stainless bucket
Please specify when ordering main unit.		Chamber temperature sensor
		Sample temperature sensor
		Pressure gauge
		External alarm output terminal
		Temperature output terminal
		Time-up output terminal